



Tactical Trauma Treatment Bandage T3

The T3 Bandage provides a multifunctional solution to a wide range of injuries in a compact package. While facilitating the treatment of conventional wounds, abdominal and large wounds, amputations, entrance/exit wounds and more, the bandage allows cleaning, packing, coverage and compression of single and multiple wounds.

This all-in-one bandage reduces the amount of equipment military personnel carry and is suitable for first aid combat emergency kits. The bandage is easily applied by untrained personnel and makes an ideal personal field dressing for every soldier.

The T3 Bandage features an fully expandable and detachable large sterile dressing pad, and elasticized cloth leader, an easy to use closure bar and an extractable gauze for cleaning/packing and secondary wound dressing.

MULTI BANDAGE FEATURES

- Multi-functional bandage for various wounds
- Expandable sterile pad – effective dressing for severe wounds, burns and amputations
- Detachable sterile pad – effective dressing for entry/exit wounds and multiple wound locations
- Folded gauze for cleaning/packing and dressing
- Compact, lightweight, waterproof, vacuum sealed packaging
- Sterile, long shelf life up to 8 years expiry
- Integral pressure applicator exerts immediate direct pressure on wound
- Built-in closure bar, no pins, no clips, no tape, no hook-and-loop, no knots
- Non-adherent pad eliminates removal pain and wound re-opening
- Compact and durable packaging

PRODUCT SPECIFICATIONS



Type. T3

Description. 4" Customizable bandage with expandable pad and extractable gauze

Elastic Bandage Size.

370cm+- 20cm length (when stretched) x 10cm width

Pad Size. 10cm width x 18cm length

Gauze Size. 75cm x 10cm width un-stretch


Bag Size. 12.5cm width x 24cm length

Part No. 105595000




T3 Additional Features

Expendable sterile pad for large wound / abdominal injury

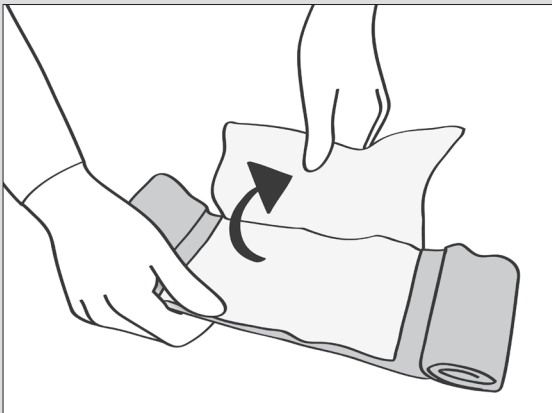

Folded gauze for cleaning /packing and dressing

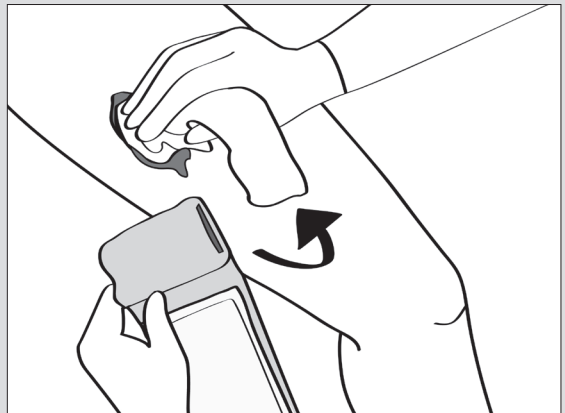

Detachable sterile pad for exit-entry wound and as additional pad for 2 separate wound location

